




Simple Operation. Precise Results.

Measuring moisture content has never been easier than with the MB90 Moisture Analyzer. The MB90's halogen heating distributes heat uniformly over the sample and speeds up the drying process, delivering fast and precise measurements. Designed for tool-free cleaning, its components can be easily removed for quick and easy maintenance. Operating the MB90 is intuitive with the icon-driven menu navigation on the touchscreen display.

Unique Features Include:

- Halogen heating increases efficiency for faster results.
- Intuitive touchscreen menu offers step-by-step guidance for any application.
- Tool-free cleaning for easy maintenance.

MB90 Moisture Analyzer

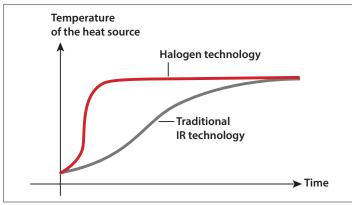

Fast Heating Time



Uniform Drying

Easy Tool-Free Cleaning

The MB90 features a patented pan handler and removable parts for convenient maintenance.

- Stainless Steel Drying Pan Wipe down quickly
- Removable Heating Chamber Components
 No tools required—simply remove glass and tray
 to clean

Exceptional Precision

The MB90 provides precise results for more confident testing.

- Fast Halogen Heating Technology
 Halogen heating dries samples in seconds and performs up
 to 40% faster than traditional infrared methods.
- Uniform Drying for Precise, Repeatable Results
 The reflective sample chamber ensures uniform heating of
 each sample to yield reliable results of 0.01% readability.
- Flexible and Powerful Features a full operating temperature range of 40–200 °C in 1 degree increments.
- Two Drying Profiles to Ensure Accuracy
 Customize heating of specific samples with temperature control for one-step drying—set temperature and time and begin!

Fast: Quickest possible time to temperature

Standard: Minimal overshoot while achieving temperatures quickly

Shut-Off Criteria for Automatic Test Completion
 Customize your own test endpoint criteria to automatically
 end the test once moisture loss is no longer detected.
 Choose from 3 criteria—Manual shut-off with pre-set time,
 Timed and Auto.



MB90 Moisture Analyzer

Testing Made Simple

The MB90 features a wide, bright touchscreen display and intuitive menu navigation for fast, accurate moisture analysis.

- Touchscreen Display
 Easy to read and understand
- Icon-Driven Menu Navigation
 Easy to follow—no training required
- 4 Mechanical Keys for Fast One-Step Operation Tare, Start/Stop, Print, On/Off
- USB and RS232 Ports
 Quick connections to a PC or flash drive
- Moisture Analysis Made Simple
 Just 3 easy steps for routine sample testing!
 - 1. Tare empty sample pan
 - 2. Add sample according to best practices
 - 3. Press start and see your results!





Multiple Ways to Connect

Versatile connectivity options to handle data collection and transfer needs—enhancing the user experience.

- USB and RS232 Ports

 Quick connections to a PC or flash drive
- Data Collection and Transfer
 You can simply save and transfer data from the moisture analyzer to a USB drive

MB90 Moisture Analyzer

Model	MB90
Capacity	90 g
Readability	0.01%/0.001 g
Repeatability (Std Dev) (g)	0.18% (3 g sample) 0.02% (10 g sample)
Moisture Range	0.01% to 100% (0.01% to 1000% for regain mode)
Output	RS232, USB Host, USB Device
Heating Element	Halogen
Temperature Range	40–200 °C
Power	240 VAC 50/60 Hz
Operating Temperature Range	10−40 °C
Display Type	109 mm, QVGA, TFT touch screen
Display Results	%moisture, %solids, %regain, time, temperature, weight, method name, drying curve and statistics
Pan Size	90 mm
Dimensions (W \times H \times D)	211 × 180 × 355 mm
Net Weight	5.23 kg
Shipping Weight	8.35 kg

Compliance

- Product Safety: IEC/EN 61010-1; IEC/EN 61010-2-010; CAN/CSA C22.2 61010-1; CAN/CSA C22.2 61010-2-010; UL 61010-1
- Electromagnetic Compatibility: IEC/EN 61326-1 Class B, Basic environments; FCC Part 15 Class A; ICES-003 Class A


Additional Features

ABS top housing, stainless steel pan support, in-use cover, storage up to 2 methods and 100 test results in total, user-guide, 14 operational languages

Accessories

Temp. Calibration Kit, MB Series	11113857
In-Use-Cover, MB90 MB120	30284478
Pan Holder, MB120&MB90	30246947
Pan, Reusable, Set(3), 7mm, MB Series 8	80252478
Pan, Reusable, Set(3), 14mm, MB Series 8	80252479
Pan, Aluminum, Set(50), MB Series 3	30585411
Filter, Glass Fiber, Set(200) 8	80850087
Cage, sample	80252477
Anti-Theft Device, MB AX EX Scout 8	80850043
Scoop, Accessory, MB9x, MB120	30284477
Printer, Impact, SF40A (EU)	30064202
Interface Kit, RS232-USB	30304101
Interface Kit, RS232-Ethernet	30304102

Outline Dimensions

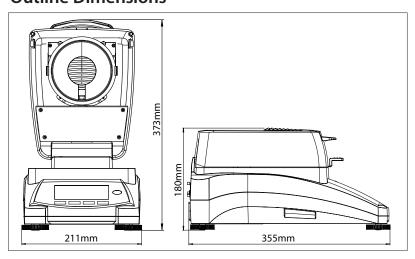

OHAUS Europe GmbH

Heuwinkelstrasse 3 8606 Nänikon Switzerland

e-mail: ssc@ohaus.com Tel: 0041 22 567 53 19 e-mail: tsc@ohaus.com Tel: 0041 22 567 53 20

www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.



80776004_G 20220411 © Copyright OHAUS Corporation