# **MEDICAL SCALES 2022**

# Wheelchair platform scale KERN MWB









# Compact wheelchair platform scale with two integrated access ramps for easy access

# Features

- Especially suitable for weighing patients in wheelchairs because of the low-profile platform which can be approached from either side
- II The large steel weighing plate also makes it ideal for weighing obese patients
- Ideal also for rapid weighing of, e.g. laundry carts, container trolleys, roller containers, etc.
- Secure and non-slip positioning with height-adjustable rubber feet
- I Level indicator to level the balance precisely
  Hold function:
- When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then read the weight
- BMI function to determine underweight/ normal weight/surplus weight
- II The scale can be easily transported using the two handles and two rollers and does not require much storage space
- Wall mount for the display device
- · Protective working cover included in delivery

#### Technical data

- Large backlit LCD display, digit height 25 mm
  Dimensions of display device
- W×D×H 200×130×60 mm
- Overall dimensions,
- W×D×H 1050×1050×100 mm • Batteries included, 6×1,5 V AA,
- operating time up to 40 h
- Cable length of display device approx. 300 cm
- Net weight approx. 55 kg

# Accessories

- Protective working cover over the display device, scope of delivery: 5 items, KERN MBC-A06S05
- Cleaning cloths, alcohol-free cleaning cloths for wipe disinfection, quick-acting, on basis of modern quaternary ammonium compounds, and also against papova viruses. Particularly gentle on materials, and therefore very well suited for the disinfection of alcohol sensitive products.
   Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540.
   Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01
- For more details see the internet



Model	Weighing range	Readability	Weighing surface	
KERN	[Max] kg	[d] kg	mm	
MWB 300K-1	300	0,1	800×800	



# PROFESSIONAL CARE

# **MEDICAL SCALES 2022**



# Pictograms



#### Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required

MEMORY

Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



To connect the balance to a printer, PC or network

Data interface RS-232:



Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.



Statistics: using the saved values, the device calculates statistical data, such as average value, standard deviation etc.

PC Software: to transfer the measurements from the device to a PC



KCP

SOFTWARE

GLP/ISO log: With date and time. Only with KERN printers

## **KERN Communication Protocol (KCP):**

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



**Piece counting:** 

Reference quantities selectable. Display can be switched from piece to weight



**Totalising level A:** 

The weights of similar items can be added together and the total can be printed out



• ગ્રેગ

TOL

Weighing units: Can be switched to e.g. nonmetric units.

Please refer to website for more details

# Weighing with tolerance range:


(Check weighing) Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model



#### Hold function:

When patients do not stand, sit or lie completely still, a stable weight is calculated using an average weight

# Your KERN specialist dealer

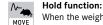

#### When the weighing conditions are unstable, a stable weight is calculated as an average



ZERO: Resets the display to "0"



Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram cf. DIN EN 60529:2000-09, IEC 60529:1989+A1:1999+A2:2013



Suspended weighing: Load support with hook on the underside of the balance



**Battery operation:** Ready for battery operation. The battery type is specified for each device



Rechargeable battery pack: Rechargeable set



Universal plug-in power supply: with universal input and optional input socket adapters for A) EU, CH; B) EU, CH, GB, USA





пмѕ

Integrated power supply unit: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, AUS or USA on request




Weighing principle: Strain gauges Electrical resistor on an elastic deforming body



Peak hold function: capturing a peak value within a measuring process



Push and Pull: the measuring device can capture tension and compression forces

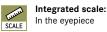

360° rotatable microscope head





MONO

Monocular Microscope: For the inspection with one eye



00

BINO

# Trinocular Microscope:

Binocular Microscope:

For the inspection with both eyes

For the inspection with both eyes and the additional option for the connection of a camera

## Abbe Condenser:

With high numerical aperture for the concentration and the focusing of light



 $\mathbf{A}$ 

ABBE

Halogen illumination: For pictures bright and rich in contrast

Ð LED

LED illumination:

Cold, energy-saving and especially long-life illumination



Fluorescence illumination for compound microscopes:

With 100W mercury lamp and filter



Fluorescence illumination for compound microscopes: With 3W LED illumination and filter




۲

∰

POLAR

Phase contrast unit: For a higher contrast

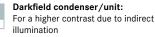

To polarise the light

Polarising unit:



Infinity system: Infinity corrected optical system



Automatic temperature compesation: For measurements between 10 °C and 30 °C



Verification possible: The time required for verification is specified in the pictogram



Package shipment: The time required for internal shipping

preparations is shown in days in the pictogram



Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram